


Multipurpose abdominal transducer

8830



- Superior image quality for all-round abdominal scanning
- Simple and convenient interventional procedures
- Comfortable design

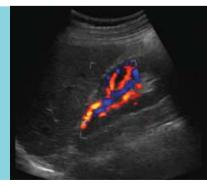

Normal kidney with flow



Normal liver

Application versatility: The 8830 is ideal for abdominal, urological, and OB/GYN scanning.

Superior image quality for all-round abdominal scanning

Scan the liver, gall bladder, pancreas and kidney with the 8830. The 8830 is the ideal all-round ultrasound tool for abdominal, urological and OB/GYN scanning. The 8830 gives you sharp images and clear details at all depths.

Simple and convenient interventional procedures

Easy-to-use needle guides. The 8830's reusable and single-use needle guides are designed specifically for performing intervention in the abdominal region, and you can continually monitor the path of the needle on the screen. The single-use needle guide is designed to guide ablative procedures.

Comfortable design

The 8830's slim design and rounded handle make it easy to position the transducer accurately, and the small foot print makes even difficult areas easy to scan. Start and stop scanning and freeze and unfreeze images with a simple press on the integrated control button.



A reusable and a single-use needle guide are available to help make intervention easier.

Specifications

Frequency range 2-6 MHz

Image field 60°

Physical data 77 x 104 x 27 mm

Contact surface (acoustic) 62.5 x 13 mm

<u>Scanning modes</u> B, M, BCFM, Doppler, Tissue Harmonic Imaging

Disinfection Immersion

Weight 180 g



World Headquarters: Mileparken 34 · DK-2730 · Herlev · Denmark Tel: (45) 4452 8100 · Fax: (45) 4452 8199 · www.bkmed.com A wholly owned subsidiary of Analogic Corporation